






Infection-free Spirometer

Accurate and Hygienic Spirometer

SPM-7

- * Convenient A4 size print-outs allow for easy chart filing
- * Grid function allows you to use low cost plain fax paper
- * Wide and Clear Color TFT LCD
- * Touch screen for easy operation
- * PRE-POST bronchodilator Comparison
- * Rotary Key Access to All Menus
- * No Cleaning
Using disposable mouthpiece patented sensing method
- Breathing Protocol
 - FVC (Forced Vital Capacity)
 - SVC (Slow Vital Capacity)
 - MVV (Maximum Voluntary Ventilation)
- Various Formulas
Morris-Polgar, ECCS-Quanjer, Kundson-ITS
- Memory for up to 200 patients
- USB Barcode Scanner & USB Memory Storage
- PC compatible through LAN
- Multi-language (English, Russian, French, German, Italian, Chinese, Spanish, Turkish, Polish, Romanian, Portuguese) available

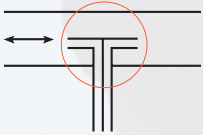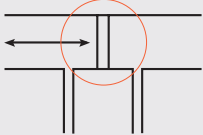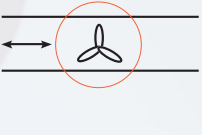

	Bernoulli's Law	Pneumotachometer	Turbine
structure			
medium	SmarTube™	screen	Propeller or turbine
accuracy	error ratio < 1% (equal air density)	error ratio < 5% (air friction)	error ratio < 20% (inertia problem)
	User-friendly Infection-free	High probability of infection because of remain saliva or bacteria.	High probability of infection because of remain saliva or bacteria and low accuracy.

SPECIFICATIONS

Dimension	47(W) × 200(H) × 34(D)mm, approx.250g
Measuring Values	FVC: FVC, FEV 1.0, FEV 1/FVC, FEF 0.2-1.2L, FEF 25-75%, FEF 75-85%, PEF, FEF 25%, FEF 50%, FEF 75%, FIVC, FEV6, PEFT, FET100% Extrapolated Volume SVC: SVC, ERV, IRV, TV, EC MVV: MVV, FB, TV
Presentation	Flow volume loop Flow time plot Measurement values table
Measuring range	Flow : 0 to ± 14 L/s Volume : 0 to ± 12L
Measuring method	Differential pressure method
Prediction equation	Morris-Polgar, ECCS-Quanjer, Knudson-ITS, Pereira, Korea CJK
Sample rate	200 samples/sec
Flow impedance	<0.2 mbar s/l at 12 l/s
Measuring accuracy	Complies with ATS (American Thoracic Society)
Environment	Ambient temperature : 15 to 40℃ (59 to 104° F) Relative humidity : 10 to 90%(non-condensing) Atmospheric pressure : 700 to 1060hPa
Standard accessory	Disposable mouthpiece 2 EA Nose clip 1 EA Adapter 1 EA Handle dock 1 EA Operation manual 1 EA Diagnosis guide 1 EA Disposable mouthpiece 1 box (100 EA)

*Specifications can be changed without prior notification.

شرکت ارکان آرا تجارت البرز

Iran : Arkan Ara Tejarat Alborz Co., Ltd.

Tel : +98-21-81060

E-mail : info@arkanara.com

www.arkanara.com

South Korea: Bionet Co., Ltd.

Tel: +82-70-7585-6418 (Int'l)/+82-2-6300-6483 (Domestic)

E-mail: sales@ebionet.com (Int'l)/ajy@ebionet.com (Domestic)

www.ebionet.com

Bionet

North America (U.S.A): Bionet America, Inc.

Tel: +1-714-734-1760

E-mail: sales@BionetUS.com

www.BionetUS.com

Latin America: Bionet Latin America, Inc.

Tel: +1-949-800-9585

E-mail: bla@ebionet.com

www.bionetesp.com